

seward

World leaders in sample preparation



# Seward Stomacher® 400 EVO



Experience the **EVO**lution of the  
Seward Stomacher®

The all new sample preparation solution  
that combines stress-free working with  
cost savings



The Seward  
Stomacher® 400 EVO

- ✓ Unique one-touch bag loading system
- ✓ Quiet and efficient blending
- ✓ Easy-clean interior
- ✓ Designed to work with our new ECO bags


## The world's leading laboratory blender has **EVO**lved

### The Seward Stomacher® 400 EVO

The all new design of the **Seward Stomacher® 400 EVO** combines the gold standard attributes of the original Stomacher® 400 with improved, easy-to-use features. The new Seward Stomacher® 400 EVO **provides maximum performance while minimising noise and consumable waste.**

#### Easy to use

- Unique 'one touch' bag loading system minimises sample bag handling, saving time and reducing the risk of sample contamination
- Simple programmable screen display enables complete control of processing, for consistent results with minimum effort

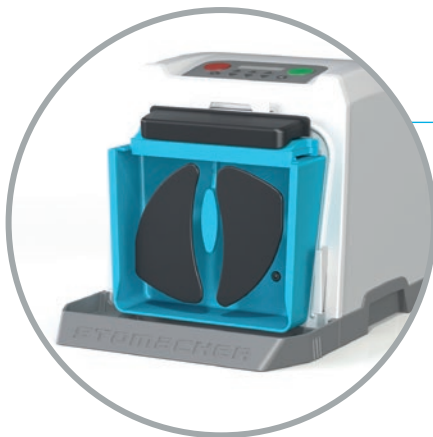

#### Easy to clean

- Front cover and strike plate can be removed to provide easy access to paddle chamber and detachable paddles for simple and thorough cleaning
- Removable drip tray allows for easy cleaning, saving time in day-to-day use

#### Quiet operation

- Double aperture seal, isolated chassis, and fully enclosed casing result in silent operation, reducing stress and disruptions to other lab work
- Rubber base with anti-vibrational technology eliminates bench movements and reduces noise



## Good sample preparation is critical for dependable results

Seward provides products that offer real value to the laboratory. Seward Stomacher® blenders have been proven over 40 years and are engineered to be reliable, strong, versatile, powerful and long-lasting. They are used to prepare over **8 million food samples in more than 150 countries** around the world, and remain the best and most trusted piece of equipment in their field.

## The gold standard in sample preparation

The Seward Stomacher® 400 EVO draws from the company's long-standing commitment to quality, reliability and performance with a sleek and durable design, supported by a 3 year warranty.

3 year warranty

## New Stomacher® 400 ECO bags

The Seward Stomacher® 400 EVO can be used with all standard bags, and you can enjoy significant cost savings by using our **new Stomacher® 400 ECO bags**.

- ✓ Save up to 25% on consumable costs
- ✓ Designed for use with the Stomacher® 400 EVO
- ✓ Significant reduction in plastic waste, ideal for ISO 14001 laboratories

Save up to **25%** on consumables\*

Stomacher® 400 **ECO** bags are designed specifically for the Stomacher® 400 EVO, and as well as saving money, **ECO** bags use up to 40% less plastic so they're more environmentally-friendly than any other blender bag.



*"Stomacher® ECO bags provide the ideal solution for streamlined labs. Less cost. Less wastage. Greater value for your lab"*

Choose from two Stomacher® 400 **ECO** bag sizes specifically suited to performing sample preparation methods for 25g Pathogen tests and 10g TVC tests.

## Stomacher® ECO bags. **ECO**nomical and **ECO**-friendly

\* Use Stomacher® 400 **ECO** bags with your Stomacher® 400 EVO and you could save up to 25% on annual consumables compared to standard Stomacher® 400 bags.

## Seward Stomacher® 400 EVO and ECO bags

The perfect solution for sample preparation



### Seward Stomacher® 400 EVO

Product Code	Description	Voltage
0400/002/EU	Stomacher 400 EVO	230V
0400/002/AJ	Stomacher 400 EVO	110V

### Stomacher® 400 ECO Bags

for Seward Stomacher® 400 EVO only

Product Code	Description	Test Volume	Dimensions (mm)	Bags per pack
BA6241	Stomacher 400 ECO bag (PATH.)	1:10 PATH Test (25g:225ml)	177 x 200	500
BA6341	Stomacher 400 ECO bag (TVC)	1:10 TVC Test (10g:90ml)	155 x 200	500

### Stomacher® 400 Standard Bags

for all Stomacher® 400 models (EVO, Circulator and Classic)

Product Code	Description	Test Volume	Dimensions (mm)	Bags per pack
BA6041	Stomacher 400 Classic Bag	400ml	177 x 305	500
BA6141	Stomacher 400 Circulator Standard Bag	400ml	177 x 305	500

\*Full range of Stomacher bags available at [www.seward.co.uk/stomacher-bags/](http://www.seward.co.uk/stomacher-bags/)

Contact Seward for further information on its range of products including the **Stomacher® 400 EVO**, and see how much you could save on your consumable costs.